



Single-Patient-Use from the beginning...

SHOP

PRODUCTS

ABOUT US

NEWS

MORE

INTL




DL

SEARCH

2022

NEODRYS®

LARGE
REFLECTIVE

 microcopy

50 LARGE
REFLECTIVE

REF DYRL

MADE IN THE USA

US 7425664

EP 1676552

Grande, Grand, Groß, Groot, Большие, 大

Ultimate Parotid Saliva Control

Lo más avanzado para el Control de la Saliva (Glándula Parótida)

Le plus avancée pour le Contrôle de la Salive (Glande Parotidée)

Ultimative Parotisspeichel-Steuerung


Lo più avanzato per il Controllo della Saliva (Ghiandola Parotide)


O Mais Avançado Controle da Saliva da Parótida


Uitmeeste speekselbeheersing

Максимальный контроль слюноотделения

優れた耳下腺唾液吸収能力

 dryifu.com

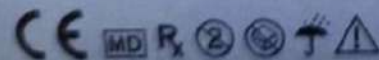

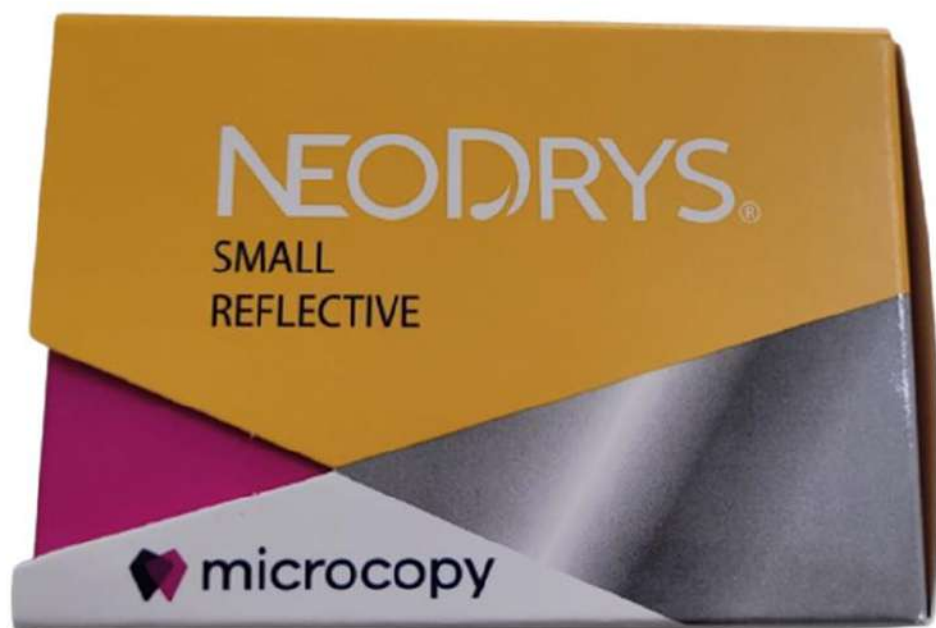

The Catalog

NEODRYS®

Suck it up buttercup

NEODRYS saliva absorbents do just that! We know how important it is to work in a dry field when doing the procedures you do all day long. Built with a super absorbent core ensures things stay nice and dry.

To serve and protect (the cheek that is)

NeoDRYS protective backing and our patented design gives you great cheek retraction. You have room to do what you need to do in close quarters, plus they protect the cheek in case you have a runaway bur situation.

Reflect on your bad choices...

And make some improvements! With a reflective backing option, NeoDRYS reflective brightens up the oral cavity and ensures you see exactly what you're working with. Out of the darkness and into the light!



ANATOMY OF ABSORPTION (HERE'S WHY WE'RE BETTER)

CORE STIFFENER

Stiffener insert provides optimum rigidity for cheek retraction.

SUPER ABSORBENT CORE

Acrylate polymer traps moisture as a gel when in contact with saliva.

GENTLE TO THE BUCCAL MUCOSA

Medical-grade poly netting adheres gently to tissue.

"QUICKNESS" COMFORT EDGE®

New technology creates a soft, flexible edge with no sharp corners while retaining product integrity.

PROTECTIVE BACKING

Choice of Original white or Reflective backing brightens the oral cavity while providing cheek protection.



WHITE BACKING

50 per pack | Small Original | Large Original



REFLECTIVE BACKING

50 per pack | Small Reflective | Large Reflective

"A must have product for day in - day out dentistry today."

Andrew Bargman DDS, Hollywood, FL

TIME IS MONEY

15 minutes of absorption saves you on both.



Instructions for Use



Color to Cheek



Point to Back



Spray to Release

Single use | microcopy.com

PRICING

Contact your dealer for pricing.

NEODRYS®

Ultimate Parotid Saliva Control

15
Minute
Parotid Control

NEODRYS Outperforms the Competition

NEODRYS advanced inner core traps moisture quickly to provide a superior dry field for up to 15 minutes. No messy leaks, even with the core fully saturated. Its moisture-proof backing protects the cheek from abrasion and bounces light into the working area. Our patented Comfort Edge® eliminates stiff abrasive material that causes patient discomfort. Absorption, protection and comfort make NEODRYS far superior to any other product on the market.

»FEATURES...




Color-Coded Sizes

NEODRYS are color coded for quick identification. Available in two sizes: Large (blue) and Small (yellow). Packaged in convenient easy-to-handle packs of 50 each.



Choice of Backing

Available in ReFlective light-enhancing, mirror-like finish or Original white backing.

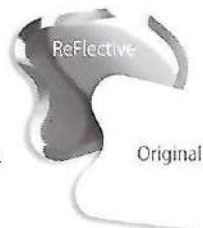

Product Codes:

NEODRYS Small Original	DYS
NEODRYS Large Original	DYL
NEODRYS Small ReFlective	DYRS
NEODRYS Large ReFlective	DYRL

 **MICROCOPY**
INNOVATIONS FOR DENTISTRY

»BENEFITS...

NEODRYS Save You Money

Stop wasting money on pads and rolls that either fall apart or get saturated quickly. When you add up the reduction in botched procedures & wasted time — NEODRYS cost less!


- ◆ No loose fibers left behind in the oral cavity.
- ◆ Saves time by protecting procedures from moisture contamination from the parotid gland for up to 15 minutes.
- ◆ Super absorbent core — far outlasts other absorbent pads.
- ◆ Complete coverage of the parotid orifice and entire buccal mucosa with a non-cutting comfortable poly barrier.
- ◆ Cheek retraction facilitates work far back in the oral cavity.
- ◆ Comfort Edge® provides maximum patient comfort.
- ◆ Stays in place yet removes easily with water spray.
- ◆ ReFlective backing provides light into the oral cavity.

No Questions Asked Guarantee

Microcopy unconditionally guarantees your total satisfaction with the quality and performance of all our products.

»APPLICATION...

Select NEODRYS size which adequately covers the buccal mucosa. Insert as shown with color side against cheek and point to back of mouth. In a few seconds, NEODRYS will begin to adhere to the tissue and stay in place. **Important:** To remove without irritating tissue, release adhesion with ample water spray to the buccal side of the NEODRYS.

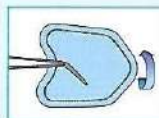

COLOR TO CHEEK



POINT TO BACK



SPRAY TO RELEASE

ORIGINAL

FLAPS[®]


Cushioned Bite Tabs



ORIGINAL




LONG

 microcopy



FLAPS[®]

Cushioned Bite Tabs
For X-Ray Digital Sensor & Film

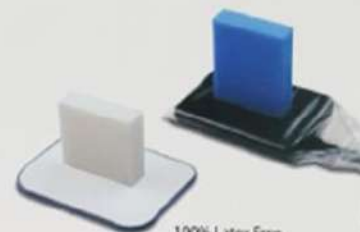

100% Latex Free
Soft Polyethylene Foam.
MADE IN THE USA

INSTRUCTIONS:



All instructions at mcifu.com.

1. Stick FLAPS on the image side of the x-ray sensor or film. Flaps will adhere to the plastic.

Required: snug-fitting sensor sheaths!

2. Position the x-ray sensor or film and instruct the patient to bite the FLAPS.
3. Ask for a raise because you've just upped your x-ray taking game.



Perfect for horizontal or vertical bitewings.



Excellent for any periapical x-ray, using either the parallel or bisecting angle technique.

FLAPS®

Simple as that!

Self-sticking foam tabs attach quickly and are easy to position.

We're on point


Cushioned tab holds film or sensor securely, allows full occlusion and provides visual reference for accurate placement and alignment.

Tooth pillows?

Soft cushioning without hard surfaces, provides the ultimate in comfort and small size reduces gagging.

Kick 'em to curb

Use-and-toss convenience—extremely cost effective tool!

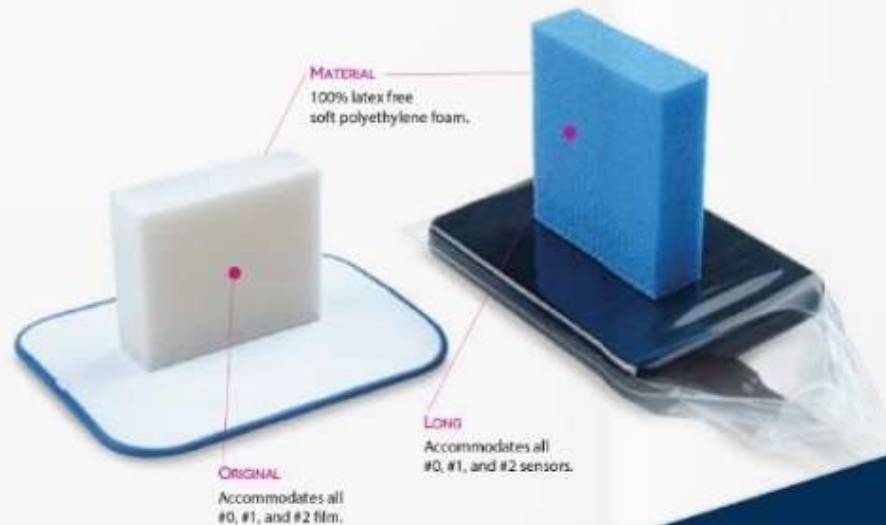

MATERIAL
100% latex free soft polyethylene foam.

ORIGINAL
Accommodates all #0, #1, and #2 film.

LONG
Accommodates all #0, #1, and #2 sensors.

PRICING
Contact your dealer for pricing.



ORIGINAL (WHITE)

500 per pack



LONG (BLUE)

375 per pack

ORIGINAL (WHITE)



Perfect for horizontal or vertical bitewings.

LONG (BLUE)



Excellent for any periapical x-ray, using either the parallel or bisecting angle technique.

Instructions for Use



Single use

Multi-use

Step 1: Stick the Flaps on the back of the film or sensor. (See diagram).

Step 2: Instruct patient to position tongue against sensor and bite onto Flaps.

Step 3: Ask for a raise because you've just upped your x-ray taking game.

Requires: Snug-fitting sensor sheaths!

"Flaps are very easy, simple and fast to use. Plus, patients don't gag like they do with other devices."

Tamara Simac, DDA

KIDS DIG US

Soft cushion keeps all patients comfortable (especially kids).

